

# REJUVE

Portable RF System for Facial Contouring & Skin Tightening Treatments.

## REJUVE FACE

### Skin Tightening & Facial Rejuvenation

As our skin matures and is continuously subjected to environmental effects, its connective tissue begins to break down, causing wrinkles and other signs of aging develop. The Multi-Polar RF energy penetrates the skin layer to selectively heat the skin tissue and the extra cellular matrix where the collagen fibers are selectively heat the skin tissue and the extra cellular matrix where the collagen fibers are embedded. The local heating causes an immediate contraction of collagen fibers, and at the same time increases the metabolism of the fibroblasts thereby accelerating the production of new collagen and elastin fibers. The contraction of collagen fibers and the regeneration of the collagen and elastin deposits, tighten the skin layers creating smoother, healthier and younger looking skin.

## FEATURES

Rejuve is painless, effective and easy to use :

- Non-invasive, relaxing treatment
- Creates a noticeable clinical effect from the very first treatment
- Significantly reduces wrinkles and fine lines
- Improves skin texture
- Long lasting anti-aging results



Rejuve is painless, effective and easy to use



## REJUVE BODY


### Body Contouring & Cellulite Treatment

- Cellulite condition: Accumulation of a large number of enlarged fat cells in the hypodermis.
- The accumulated fat presses the skin creating the cottage cheese skin effect.
- Multipolar RF technology selectively increases the temperature of fat cells doing so simultaneously in deep and superficial layers of fat.
- This increases metabolism and release of liquid fat from the fat cells.
- Normal cells can support 60 , but other cells including tumor cells or fat cells burn or come to apoptosis on 55
- Fat starts to be melted at 41 , and contains numerous venous-lymphatic vessels which allow a strong electromagnetic field. Moreover, high reflection in the frontier of fat layer and muscle layer leads easily the diathermy in the fat.

We can share our success stories with **Rejuve**

## SPECIFICATION

- Model Name : Rejuve
- Power Input : 100/230V, 50/60Hz
- Max. Power Consumption : 30W
- Dimensions : 32(W) 25(D) 15(H)cm
- Weight : 5kg
- Output Frequency : 1MHz
- Output Mode : Three-Polar
- Treatment Time : 60min max.
- Safety Device : Zero Start, Safety Switch
- Intensity Control : PWM & Swing

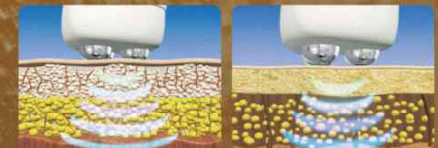

DAEYANG MEDICAL CO., LTD.

#1656-1, Donghwari, Munmakeup, Wonju-si, Gangwon-do, 220-801, KOREA.

• Tel : +82-33-747-4466 • Fax : +82-33-747-1133 • e-mail : daeyang@daeyangmed.com • website : www.daeyangmed.com